


a helpful
guide to
our product
offerings.



PENTRON

Dentistry Open Wide™

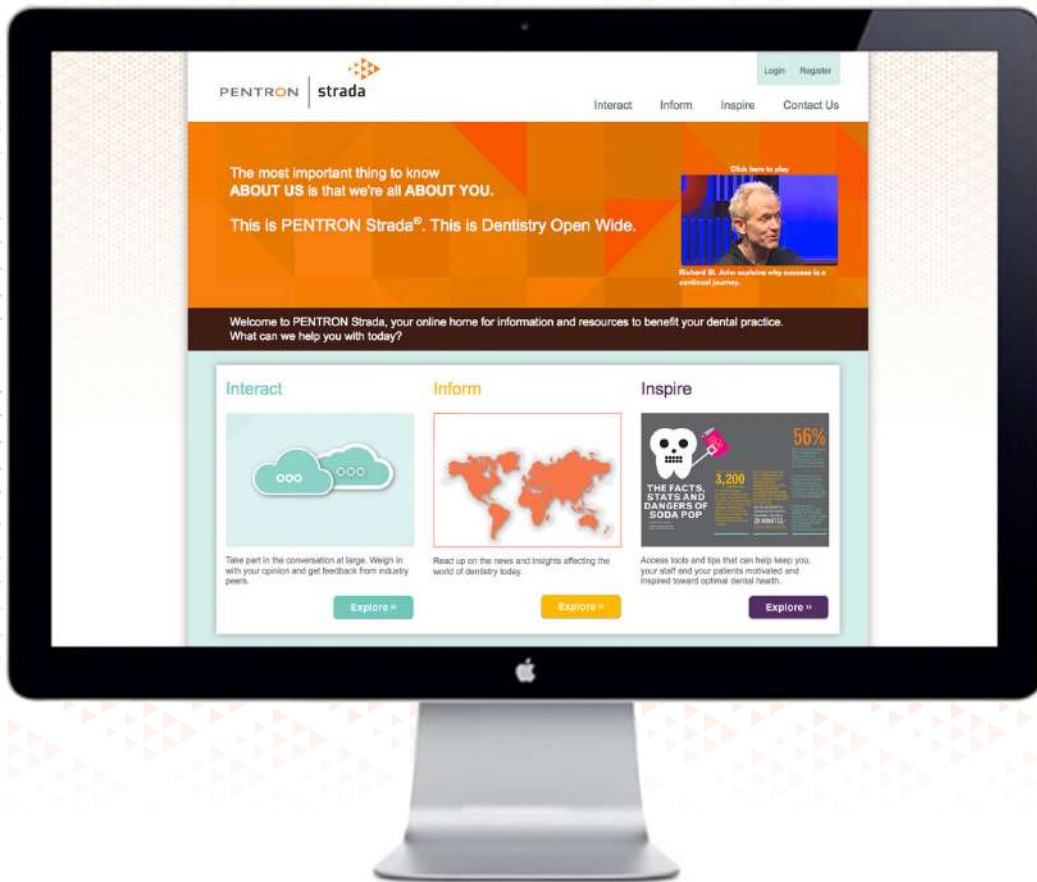

So What's the Story Behind PENTRON Strada™?

As a manufacturer of dental products, PENTRON has a long-standing reputation for delivering recognizable quality and value in all they do. Like the dental professionals they serve, PENTRON understands that it takes more than just great tools to achieve optimal dental health. So, to demonstrate this belief and remain true to their heritage, PENTRON wanted to create deeper, more comprehensive relationships with the dental professionals they serve; a relationship built less on the extensive PENTRON product line and more on creating solutions that can help dentists build successful practices.

The name for this ambitious new venture is PENTRON Strada, which was derived from the concept of “a new path, road, or means to a new relationship” between manufacturer, customers, and end-users. It lives as an online platform at PentronStrada.com, where it can create meaningful dialogue between dental professionals and PENTRON, fill information gaps within the practice, and provide problem-solving strategies. Databases within PentronStrada.com are fed by incoming news feeds and user interaction, and can grow in depth and breadth as time goes on. In other words, the more the site is used, the more it can offer its users.

As always, PENTRON continues to sell dental restorative materials through the parent site: Pentron.com. Meanwhile, PentronStrada.com rounds out the story—positively impacting dental practices with information and tools that have nothing to do with the PENTRON product line. PentronStrada.com is proof that dental professionals can always get something valuable from PENTRON—even without making a purchase.

TABLE OF CONTENTS

COMPOSITES

Fusio™ Self-Adhesive Flowable Composite	p. 3
Flow-It™ ALC Flowable Composite	p. 4
Alert™ Condensable Composite	p. 5
Artiste™ Nano Composite.....	p. 6
Simile™ Universal Composite	p. 7

POST & CORE

FibreKleer™ 4x Fiber Post	p. 8
FibreKor™ Post	p. 9
Build-It™ FR Core Buildup Material	p. 10
Build-It™ Light Cure Core Material.....	p. 11

TEMPORARY

TempSpan™ Temporary C&B Material.....	p. 12
TempSpan™ Clear Matrix Material.....	p. 12
TempSpan™ Transparent Temporary Cement.....	p. 13
Zone™ Non-Eugenol Temporary Cement & ZONEfree Translucent Temporary Cement.....	p. 14-15
Embonte Temporary Cement.....	p. 16

IMPRESSION

Correct Plus™ Impression Material	p. 17
Correct Plus™ Fast Set Impression Material	p. 18
Correct Plus™ & Correct Quick™ Bite Registration Material & Correct VPS™ Putty.....	p. 19
Image Dust Free Alginate.....	p. 20
VanR Reversible Hydrocolloid Impression Material	p. 21-27

CEMENTS

Mojo™ Veneer Cement.....	p. 28
Breeze™ Self-Adhesive Resin Cement.....	p. 29
Cement-It™ Universal C&B Resin Cement.....	p. 30

BONDING

Bond-1™ Primer/Adhesive	p. 31
Bond-1™ SF Solvent-Free SE Adhesive	p. 32
Bond-It™ Bonding System	p. 33

ACCESSORIES

Etching Gel, Adhesives Accessories, Impression Accessories	p. 34
Impression Trays	p. 35
Tip Refills.....	p. 36-37
Dispenser Guns.....	p. 38

COMPOSITES

Fusio™ Self-Adhesive Flowable Composite

Fusio Self-Adhesive Flowable Composite represents the natural progression of composite resin technology that clinicians have been waiting for. *Fusio* combines the benefits of adhesive and restorative technology into one product bringing restorative techniques to new unprecedented heights. Simply syringe into the preparation, agitate, and light-cure. Immediately upon application *Fusio* conditions, bonds, and seals without the need for an etchant or an adhesive. With such unique applications as dentin replacement and self-adhesive base liner, *Fusio* will change the way you look at restorative dentistry. Best of all *Fusio* costs less than premium flowable composites and their required separate adhesives.

Features & Benefits:

- **Self-adhesive** – bonds to dentin and enamel without a separate bonding agent
- **Radiopaque** – easily stands out on an x-ray for quick and easy diagnosis
- **Seals dentin** – offers the potential for reduced sensitivity
- **Multiple Indications** – pit & fissure sealant, block out undercuts, base liner, small Class I, Class III, and Class V



Fusio Self-Adhesive Flowable Composite	Includes 2 - 1ml/1.8gm Syringes, 10-19 Gauge Tips, Instructions
---	--

SHADE	REF/ PART NO.
A1	N21SA
A2	N21SB
A3	N21SC
B1	N21SD

Tip Refills	
--------------------	--

DESCRIPTION	REF/ PART NO.
50 - 19 Gauge Tips	N28NT

COMPOSITES



Flow-It™ ALC Flowable Composite

Flow-It ALC (Accelerated Light Cure) is a versatile flowable composite with flexural and tensile strengths comparable to the best hybrid composites. Excellent handling capabilities and just the right viscosity make it easy to use in an extensive range of applications. *Flow-It ALC* cures in only 3 seconds with high-speed curing lights and in only 10 seconds with standard halogen curing lights.

Features & Benefits:

- **Reliable handling** – material flows well without slumping to make placement easy
- **Versatile** – effectively used in a wide variety of procedures such as liners for Class II, Class V restorations, repair temporaries, and small Class I & III restorations
- **Shade selection** – complete esthetic selection of 24 shades

REF N11WA *Flow-It ALC* Intro Kit

6 - 1ml/1.5gm Syringes *Flow-It ALC* shades:

A1, A2, A3, B2, C3, Universal Opaque

- Instructions, 50 - 20 Gauge Tips

Flow-It ALC Flowable Composite Syringe Refills

Includes 1 - 1ml/1.5gm Syringe, Instructions, 5-20 Gauge Tips

SHADE	REF/ PART NO.	SHADE	REF/ PART NO.	SHADE	REF/ PART NO.
A0	N11XA	C1	N11N	Univ. Opaque	N11I
A1	N11A	C2	N11F	Incisal	N11XC
A2	N11B	C3	N11G	Pedo XLgt.	N11V
A3	N11C	C4	N11O	Pedo Lgt.	N11U
A3.5	N11J	Clear	N11H	White Opaque	N11XD
A4	N11K	D2	N11P		
B0	N11XB	D3	N11Q		
B1	N11L	D4	N11R		
B2	N11D				
B3	N11E				
B4	N11M				

Flow-It ALC Flowable Composite Value Packs

Includes 6 - 1ml/1.5gm Syringes of one shade, Instructions, 50-20 Gauge Tips

SHADE	REF/ PART NO.	SHADE	REF/ PART NO.	SHADE	REF/ PART NO.
A0	N11VXA	C1	N11VN	Univ. Opaque	N11VI
A1	N11VA	C2	N11VF	Incisal	N11VXC
A2	N11VB	C3	N11VG	Pedo XLgt.	N11VV
A3	N11VC	C4	N11VO	Pedo Lgt.	N11VU
A3.5	N11VJ	Clear	N11VH	White Opaque	N11VXD
A4	N11VK	D2	N11VP		
B0	N11VXB	D3	N11VQ		
B1	N11VL	D4	N11VR		
B2	N11VD				
B3	N11VE				
B4	N11VM				

Tip Refills

DESCRIPTION	REF/ PART NO.
50 - 20 Gauge Tips	T26A
100 - 20 Gauge Tips	T26C

COMPOSITES

Alert™ Condensable Composite

Alert Condensable Composite is a posterior restorative material designed to be a viable amalgam alternative. In addition to over 10 years of independently validated* clinical success, what sets *Alert* apart from other composites is its proprietary chopped glass fiber filler that provides amalgam-like handling and exceptional physical properties.

Features & Benefits:

- **Efficient** – low polymerization shrinkage of 1.9%
- **Long-lasting** – minimal wear rate of only 3.5µm per year
- **Easy-to-use** – handles and packs similar to amalgam
- **Accommodating** – available in the 5 most popular posterior shades: A2, A3, A3.5, B1 & C2



Alert Condensable Composite Syringe Refills	Includes 1 - 4gm Syringe, Instructions
---	--

SHADE	REF/ PART NO.	SHADE	REF/ PART NO.
A2	N15SA	B1	N15SD
A3	N15SB	C2	N15SE
A3.5	N15SC		

Alert Condensable Composite Single-Use Refill	Includes 10 - 0.25gm Single Doses, Instructions
---	---

SHADE	REF/ PART NO.	SHADE	REF/ PART NO.
A2	N15HA	B1	N15HD
A3	N15HB	C2	N15HE
A3.5	N15HC		

Accessories	
-------------	--

DESCRIPTION	REF/ PART NO.
Single Dose Dispenser Gun	N75A

COMPOSITES

Artiste™ Nano Composite

Artiste Nano Composite is a light cured, radiopaque, nano-engineered restorative material specially formulated to adapt to all anterior and posterior cavity preparations. Easy to handle, Artiste combines esthetics and strength for exceptional polishability and performance – all in one universal composite.



Features and Benefits:

- **Ideal handling** – non-slumping, non-sticky consistency for effortless sculpting and shaping
- **Nano-engineered chemistry** – the perfect combination of esthetics and strength for exceptional polishability and performance
- **Intuitive shading system** – follows the Vita® shade guide for easy reference and includes bleach & A opaque shades
- **Convenient delivery system** – choice of single dose or syringe

Artiste Nano Composite Syringe Refills

Includes 1 - 4gm Syringe, Instructions

SHADE	REF/ PART NO.	SHADE	REF/ PART NO.	SHADE	REF/ PART NO.
A1	N36AA	B2	N36AG	Bleach	N36AQ
A2	N36AB	B3	N36AH	A Opaque	N36AOA
A3	N36AC	C1	N36AJ		
A3.5	N36AD	C2	N36AK		
B1	N36AF	C3	N36AL		

Artiste Nano Composite Single-Use Refills

Includes 10 - 0.25gm Single Doses of One Shade, Instructions

SHADE	REF/ PART NO.	SHADE	REF/ PART NO.	SHADE	REF/ PART NO.
A1	N36BA	B2	N36BG	Bleach	N36BQ
A2	N36BB	B3	N36BH	A Opaque	N36BOA
A3	N36BC	C1	N36BJ		
A3.5	N36BD	C2	N36BK		
B1	N36BF	C3	N36BL		

Accessories

DESCRIPTION	REF/ PART NO.
Single-Dose Dispenser Gun	N75A

COMPOSITES



Simile™ Universal Composite

Simile Universal Composite combines strength and esthetics for both anterior and posterior restorations. Its nano-hybrid formula has been independently proven* to deliver superior surface smoothness for easy polishability, superior polish retention, exceptional wear resistance, and outstanding physical properties.

Features & Benefits:

- **Sculptable** – handling that will not slump and will not stick to your instrument
- **Polishability** – nano-hybrid fillers produce a simple polish that lasts over time
- **Optimum shade selection** – available in 20 shades to meet all your esthetic needs

REF N56A **Simile 10 Shade Syringe Kit**

10 - 4gm Syringes

Shades: A1, A2, A3, A3.5, B1, B2, B3, C2, Opaque Universal, Clear Incisal

- Instructions

REF N56C **Simile 10 Shade Single Dose Kit**

100 - 0.25gm Single Doses

Shades: A1, A2, A3, A3.5, B1, B2, B3, C2, Opaque Universal, Clear Incisal

- Instructions

Simile Universal Composite Syringe Refills

Includes 1 - 4gm Syringe, Instructions

SHADE	REF/ PART NO.	SHADE	REF/ PART NO.
A0	N56AT	C1	N56AL
A1	N56AA	C2	N56AH
A2	N56AB	C3	N56AK
A3	N56AC	C4	N56AP
A3.5	N56AD	D2	N56AQ
A4	N56AN	D3	N56AR
B0	N56AU	D4	N56AS
B1	N56AE	Opaque Universal	N56AI
B2	N56AF	Clear Incisal	N56AJ
B3	N56AG		
B4	N56AO		

Simile Universal Composite Single-Use Refills

Includes 15 - 0.25gm Single Doses of one shade, Instructions

SHADE	REF/ PART NO.	SHADE	REF/ PART NO.
A0	N56BS	C1	N56BH
A1	N56BA	C2	N56BI
A2	N56BB	C3	N56BJ
A3	N56BC	C4	N56BO
A3.5	N56BD	D2	N56BP
A4	N56BM	D3	N56BQ
B0	N56BT	D4	N56BR
B1	N56BE	Opaque Universal	N56BK
B2	N56BF	Clear Incisal	N56BL
B3	N56BG		
B4	N56BN		

POST AND CORE



FibreKleer™ 4x Fiber Post

FibreKleer 4x Fiber Posts have outstanding radiopacity of up to 400% of aluminum — to help ensure that they are clearly identified on a radiograph at placement and follow-up appointments. Translucent glass fibers are bundled in a proprietary clear resin matrix to provide reliable flexural strength and help preserve the desired, natural lifelike esthetics. A dentin-like flexural modulus results in a post that flexes with the tooth and transmits stress more evenly along the remaining root structure to help minimize the instance of root fractures, compared to metal posts. For convenience, the Post System Kit packaging features a bur block to securely hold posts and drills while also improving accessibility.

Features and Benefits:

- **Radiopacity up to 400% of Al** – the post is clearly identified on a radiograph
- **Multiple sizes and body styles** – enables you to fit the post to the tooth—not the tooth to the post
- **Mimics properties of the root** – ideal flexural modulus properties help minimize root fractures
- **Excellent light translucency** – glass fibers help to preserve the natural lifelike esthetics



REF N83A *FibreKleer 4x Parallel Post System Kit*

15 - *FibreKleer 4x* Parallel Posts: 5 each in sizes 1.00mm, 1.25mm, and 1.50mm
3 Drills - 1 each in sizes 1.00mm, 1.25mm, and 1.50mm



REF N83B *FibreKleer 4x Tapered Post System Kit*

15 - *FibreKleer 4x* Tapered Posts: 5 each in sizes 1.25mm, 1.375mm, and 1.50mm
3 Drills - 1 each in sizes 1.25mm, 1.375mm, and 1.50mm



REF N83C *FibreKleer 4x Original Post System Kit*

15 - *FibreKleer 4x* Original Posts: 5 each in sizes *1.00mm, 1.25mm, and 1.50mm
3 Drills - 1 each in sizes 1.00mm, 1.25mm, and 1.50mm

FibreKleer 4x Post 10pk Refills Includes Instructions

PARALLEL POST	REF/ PART NO.	TAPERED POST	REF/ PART NO.	ORIGINAL POST	REF/ PART NO.
1.00mm Post (red)	N83AA	1.25mm Post (black)	N83BA	*1.00mm Post (red)	N83CA
1.25mm Post (yellow)	N83AB	1.375mm Post (purple)	N83BB	1.25mm Post (yellow)	N83CB
1.50mm Post (blue)	N83AC	1.50mm Post (green)	N83BC	1.50mm Post (blue)	N83CC

FibreKleer 4x Post 30pk Refills Includes Instructions

PARALLEL POST	REF/ PART NO.	TAPERED POST	REF/ PART NO.	ORIGINAL POST	REF/ PART NO.
1.00mm Post (red)	N83AD	1.25mm Post (black)	N83BD	*1.00mm Post (red)	N83CD
1.25mm Post (yellow)	N83AE	1.375mm Post (purple)	N83BE	1.25mm Post (yellow)	N83CE
1.50mm Post (blue)	N83AF	1.50mm Post (green)	N83BF	1.50mm Post (blue)	N83CF

FibreKleer 4x and *FibreKor* Fiber Post Drill Refills

TAPERED DRILL	REF/ PART NO.	PARALLEL / ORIGINAL DRILL	REF/ PART NO.	PARALLEL / ORIGINAL DRILL	REF/ PART NO.
1.25mm Drill (black)	N81CA	1.00mm (red)	N18DA	1.125mm (purple)	N18DD
1.375mm Drill (purple)	N81CB	1.25mm (yellow)	N18DB	1.375mm (green)	N18DE
1.50mm Drill (green)	N81CC	1.50mm (blue)	N18DC		

*Due to size, the 1.00 mm post is not serrated

POST AND CORE



FibreKor™ Fiber Post

The *FibreKor* Fiber Post system provides five post sizes to accommodate a range of canal sizes. Designed from uniquely treated glass fibers bundled in a strong resin matrix, the *FibreKor* Fiber Post system generates an integrated bond among tooth structure, bonding agent, resin cement, and composite core material.

Features & Benefits:

- **Effective** – unique design allows for excellent adherence of resin cement
- **Flexible** – easy shaping and cutting to size for any canal
- **Variety** – three sizes in the kit: *1.00mm, 1.25mm, and 1.50mm diameters with refills available in five sizes

REF N18

FibreKor Fiber Post Kit

30 - *FibreKor* Posts:

10 each of sizes *1.00mm, 1.25mm, and 1.50mm

3 Drills - 1 each of sizes 1.00mm, 1.25mm, and 1.50mm

- Instructions, Accessories

FibreKor Fiber Post 10pk Refills

Includes Instructions

SIZE	REF/PART NO.	SIZE	REF/ PART NO.	SIZE	REF/PART NO.
*1.00mm Post (red)	N18BA	1.25mm Post (yellow)	N18BB	1.50mm Post (blue)	N18BC
1.125mm Post (purple)	N18BD	1.375mm Post (green)	N18BE		

FibreKor Fiber Post 30pk Refills

Includes Instructions

SIZE	REF/PART NO.	SIZE	REF/ PART NO.	SIZE	REF/PART NO.
*1.00mm Post (red)	N18AB	1.25mm Post (yellow)	N18AD	1.50mm Post (blue)	N18AF
1.125mm Post (purple)	N18AC	1.375mm Post (green)	N18AE		

FibreKor Fiber Post Refill Combo Pack

Includes Instructions

DESCRIPTION	REF/PART NO.
10 each of sizes *1.00mm, 1.25mm, 1.50mm (red, yellow, blue)	N18A

FibreKor and *FibreKleer* 4x Fiber Post Drill Refills

SIZE	REF/PART NO.	SIZE	REF/PART NO.	SIZE	REF/PART NO.
1.00mm (red)	N18DA	1.25mm (yellow)	N18DB	1.50mm (blue)	N18DC
1.125mm (purple)	N18DD	1.375mm (green)	N18DE		

POST AND CORE

Build-It™ FR Core Material

Build-It FR Fiber Reinforced Core Material is the perfect material to use with *FiberKleer 4x* and *FiberKor* Fiber Posts. Formulated to flow and stack without slumping, the material's unique viscosity enables the clinician to dispense the material and form the core build-up with ease. Specially treated glass fibers provide a minimum compressive strength of 300MPa for a strong and reliable foundation for your indirect restorations. Designed to offer optimal versatility and performance, *Build-It FR* Core Material self cures in minutes or light cures in seconds to a hard consistency. Choose from shades A2, A3, Opaceous White, Gold, and Blue in either 25ml auto-mix cartridges or 4ml Mini Mix™ syringes.

Features & Benefits

- **Easy handling** – flows and stacks where you want it without slumping
- **Versatile** – designed to offer optimal versatility and performance
- **Proven performance** – strong and reliable foundations are clinically tested
- **Award-winning** – efficient results have won numerous awards

REF N32 *Build-It FR* Core Material Kit

- 1 - 25ml/48gm Auto-mix Cartridge A2
- 1 - *Build-It FR* 1:1 25ml Dispenser Gun
- 20 - Mixing Tips
- 20 - Intra-oral Tips
- Instructions

REF N32F *Build-It FR* Mini-Mix Core Material Kit

- 5 - 4ml Syringes
- Shades: A2, A3, Blue, Gold, and Opaceous White
- 25 - Mixing Tips
- 25 - Intra-oral Tips
- Instructions



Build-It FR Core Material - 25ml Refills

Includes 1 - 25ml/48gm Auto-mix Cartridge, Instructions, 20 Mixing Tips, 20 Intra-oral Tips

SHADE	REF/PART NO.	SHADE	REF/ PART NO.	SHADE	REF/PART NO.
A2	N32AA	Blue	N32AC	Opaceous White	N32AD
A3	N32AE	Gold	N32AB		

Build-It FR Core Material - Mini-Mix Refills

Includes 4 - 4ml Mini-Mix Syringes, Instructions, 20 Mixing Tips, and 20 Intra-oral Tips

SHADE	REF/PART NO.	SHADE	REF/ PART NO.	SHADE	REF/PART NO.
A2	N32FA	Blue	N32FD	Opaceous White	N32FE
A3	N32FB	Gold	N32FC		

Build-It FR Core Accessories

25ML CARTRIDGE ACCESSORIES	REF/PART NO.	MINI-MIX ACCESSORIES	REF/PART NO.
48 - Intra-oral Tips - yellow	N32D	48 - Mixing Tips/48 - Intra-oral Tips	N32GA
48 - Wastesaver Tips - blue	N32E	100 - Mixing Tips/100 - Intra-oral Tips	N32GB
1 - 25ml 1:1 Dispenser Gun	N32B	1 - 1:1 Auto-Mix Syringe Dispenser Gun	N75S

POST AND CORE

Build-It™ Light Cure Core Material

Build-It Light Cure Core Material is specifically designed for clinicians who favor the on-command cure afforded by light cure core materials. The high translucency of *Build-It Light Cure* enables a 10mm depth of cure without the need to build up incremental layers – with as little as 20 seconds of curing time per surface. Furthermore, once *Build-It Light Cure* is polymerized, it behaves just like original *Build-It FR* Core Material, with outstanding physical properties and the ability to prep like dentin without ditching. A proprietary BisGMA-free resin promotes non-sticky, sculptable handling that allows for quick and easy adaption to the tooth and post.

Features & Benefits:

- **10mm depth of cure** – streamlines placement and shortens procedure time
- **Easy-to-use** – handles and packs like a composite with a deeper curing option for core build-ups
- **Compatible with any bonding agent** – enables you to use your preferred bonding agent
- **Convenient delivery** – available in syringe and single dose delivery options

REF N32T **Build-It Light Cure Core Material Syringe Kit**

- 2 - 2.2ml/5gm Syringes
• Instructions, Accessories

REF N32U **Build-It Light Cure Core Material Single-Dose Kit**

- 30 - 0.25gm Single Doses
• Instructions, Accessories



TEMPORARY MATERIALS



TempSpan™ Temporary Crown & Bridge Material

TempSpan Temporary C&B Material is a dual cure composite specifically formulated to produce highly esthetic and durable provisional restorations in the least amount of time possible. TempSpan Temporary C&B Material's unique dual-cure formula allows you to light cure the provisional both intraorally and extraorally, minimizing both chair time and patient discomfort.

Features & Benefits:

- **Outstanding physical properties** – ensures durable provisional restorations
- **Dual cure formulation** – saves time and minimizes patient discomfort
- **Wide range of shades** – seven shades ensure you never have to compromise esthetics
- **Convenient** – works with standard 50ml impression dispensing gun

TempSpan C&B Material

Includes 1 - 50ml/80gm Auto-mix Cartridge, Instructions, 15 Mixing Tips

SHADE	REF/PART NO.	SHADE	REF/ PART NO.	SHADE	REF/PART NO.
A0	N69AA	A3	N69AD	B1	N69AF
A1	N69AB	A3.5	N69AE	C2	N69AG
A2	N69AC				

TempSpan C&B Accessories

DESCRIPTION	REF/PART NO.	DESCRIPTION	REF/PART NO.
30 - C&B Material Mixing Tips	N69AH	100 - C&B Material Mixing Tips	N69AI



TempSpan™ Clear Matrix Material

TempSpan Clear Matrix Material is a clear pre-operative vinyl polysiloxane impression material that is ideal when creating templates for direct composite and provisional restorations. When used in conjunction with our clear impression trays the translucency of this material permits intra-oral polymerization of any light cure composite or TempSpan Dual Cure Temporary Crown & Bridge Material. With an oral set time of 2 minutes, TempSpan Clear Matrix Material is a significantly faster alternative to vacuum-formed splints, alginate, and hand-mixed putties.

Features & Benefits:

- **Clear vinyl polysiloxane formula** – enables intra-oral polymerization minimizing chair time
- **Time-saver** – working time of 1 minute and an oral set time of 2 minutes shorten provisional procedures
- **Auto-mix delivery system** – eliminates messy, time-consuming hand mixing

REF Q34CB TempSpan Clear Matrix Material

- 4 - 50ml/52gm Cartridges of TempSpan Clear Matrix Material
- 16 - Mixing Tips
- Instructions

TempSpan Clear Matrix Accessories

DESCRIPTION	REF/PART NO.
100 - Clear Matrix Material Mixing Tips	Q19BB

TEMPORARY MATERIALS



TempSpan™ Transparent Temporary Cement

TempSpan Transparent Temporary Cement offers the perfect solution to conveniently cement your temporary restorations at a price that makes sense. Featuring eugenol-free resin chemistry, *TempSpan* Transparent Temporary Cement offers optimum bond strengths to prevent leakage and sensitivity, but can be easily removed when needed without damaging the preparation. Auto-mix syringe delivery provides hassle-free, precise placement directly into the restoration with minimal waste.

Features & Benefits:

- **Transparent** – ideal for the highly esthetic anterior zone
- **Ideal retention** – securely retains provisional restorations and facilitates easy removal
- **Non-eugenol** – ensures compatibility with resin-based permanent cements

REF N69CB *TempSpan* Transparent Temporary Cement

- 1 - 6gm Auto-mix Syringe of *TempSpan* Transparent Temporary Cement
- 12 - Mixing Tips
- Instructions

TEMPORARY MATERIALS



Zone™ – Non-Eugenol Temporary Cement

Zone is a zinc oxide non-eugenol temporary cement. Zone is available in the standard dentin color or in an A1 shade.

Features & Benefits:

Zone Automix Syringe (4g)

- Precise mixing of base and catalyst to maximize material characteristics
- Dispensed directly into any temporary
- Immediately Seated
- No mixing, no waiting, no hassle

Zone Dual Barrel Syringe (15g)

- Eliminates tubes and measuring

Zone Single Unit Dose (25 x 0.75g)

- Individual pouches for ultimate convenience

Zone Automix Syringe

Includes 1 - 4g Syringe, 8 Mixing Tips

DESCRIPTION	REF/PART NO.
ZONE Automix Syringe	27047
ZONE Automix Syringe A1	27046

Zone Dual Barrel Syringe

Includes 1 - 15g Syringe

DESCRIPTION	REF/PART NO.
ZONE Syringe	27043
ZONE Syringe A1	27045

Zone Single-Use

Includes 1 - Box of 25 Pouches

DESCRIPTION	REF/PART NO.
ZONE Single-Use	27039

Zone Kit

Includes - Kit of 1 tube base, 1 tube catalyst, and 1 non-slip mixing pad

DESCRIPTION	REF/PART NO.
Zone Kit	27040

TEMPORARY MATERIALS



ZONEfree – Translucent Temporary Cement

ZONEfree uses prismatic color blending nanofillers to blend with any adjacent surface making it the perfect choice for all your temporary cementing needs.

Features & Benefits:

- **Free of Triclosan**
- **Prismatic Nanofillers** – Unique color blending
- **For Color Blending**
- **Ideal** – Retention and strength
- **Effortless clean-up**
- **Zinc oxide non-eugenol** – non-resin formula avoids any resin-on-resin issues

ZONEfree Automix Syringe

Includes 1 - 3.2g Syringe, 8 Mixing Tips

DESCRIPTION

ZONEfree Automix Syringe

REF/PART NO.

27041

ZONEfree Dual Barrel Syringe

Includes 1 - 11.3g Syringe

DESCRIPTION

ZONEfree Syringe

REF/PART NO.

27042

ZONEfree Single-Use

Includes 1 - Box of 25 packets

DESCRIPTION

ZONEfree Single-Use

REF/PART NO.

27029

TEMPORARY MATERIALS



Embonte – Temporary Cement

Embonte Zinc Oxide Eugenol Temporary Cement is now available in a convenient automix delivery system.

Features & Benefits:

- **Mix Time** – 30-seconds
- **Fast Action** – Working time of 2 minutes at 72°F
- **Trim the Excess** – 2 minutes after seating in the mouth
- **Final Set** – 2.5 minutes after seating

Embonte Automix Syringe

Includes 1 - 4g Syringe, 8 Mixing Tips

DESCRIPTION

REF/PART NO.

Embonte Automix Syringe

27203

Embonte Dual Barrel Syringe

DESCRIPTION

REF/PART NO.

15g Syringe - Dual Barreled

27200

IMPRESSION



IMAGE Dust-Free Alginate

IMAGE Dust-Free Alginate is a high-quality alginate at a value price. *IMAGE* is a green colored, mint-scented alginate that delivers a smooth, surface with detail. Once you try *IMAGE*, you won't go back to your old alginate. *IMAGE* Dust-Free Regular Set Alginate offers a 30-second mix time, 2:20 minute working time and a 3:30 set, while *Image* Dust-Free Fast Set Alginate offers a 30-second mix time, 1:45 minute working time and a 2:20 set time.

Features & Benefits:

- **Eliminate the Mess** – Dust free
- **Bright Green Color** – Makes all details highly visible
- **Fresh Mint Odor** – Neutral scent for all

IMAGE Dust-Free Alginate	Includes 1lb. pouch, canister, scoop and water vial
DESCRIPTION	REF/PART NO.
Starter Pak, Regular Set	27426
Starter Pak, Fast Set	27416

IMAGE Dust-Free Alginate	Includes 10, 1lb bags 1 canister, 1 powder scoop, and 1 water vial
DESCRIPTION	REF/PART NO.
Value Pak, Regular Set	27427
Value Pak, Fast Set	27417
Institutional Pak, Fast Set	Includes 30, 1lb bags 27460
Institutional Pak, Regular Set	Includes 30, 1lb bags 27461

IMPRESSION



Correct Plus™ Fast Set Impression Material

Correct Plus Fast Set Impression Material is ideal when a faster-setting impression material is desired for 1-3 unit restorations. *Correct Plus Fast Set* exhibits outstanding tear strength and dimensional stability to help deliver remarkably detailed and precise impressions. See the Impression Material Comparison Chart on page 21 for a complete list of work and set times.

Features & Benefits:

- **Fast work and set times** – With min work time of 45 seconds and a total set time of 2:30 minutes, it's ideal for 1-3 unit impressions
- **High tear strength** – outstanding tear strength ensures reliable impressions
- **Versatility** – three cartridge viscosities available (Light, Medium, Heavy)

• 4 - 50ml Cartridge Pack, Instructions, 16 Mixing Tips*

• 32 - 50ml Cartridge Bulk Pack, Instructions, 132 Mixing Tips*

Correct Plus Fast Set Impression Cartridges

DESCRIPTION	4 CARTRIDGES	32 CARTRIDGES
Thick'n Thin™ Light Body	REF/PART NO. Q36C	REF/PART NO. Q36CC
Universal™ Body	REF/PART NO. Q36B	REF/PART NO. Q36BB
Auto-Mix Putty	REF/PART NO. Q36A	REF/PART NO. Q36AA

IMPRESSION

Correct Plus™ Impression Material

Correct Plus Impression Material is a vinyl polysiloxane impression material offered in a wide range of viscosities to accommodate all impression techniques. Correct Plus Impression Material is specially formulated to provide exceptional tear strength with outstanding hydrophilicity for more accurate impressions in the presence of fluids. See the Impression Material Comparison Chart on page 21 for a complete list of work and set times.



Features & Benefits:

- **Reliable** – exceptional physical properties deliver a good impression in one attempt
- **Hydrophilic** – works efficiently in the presence of moisture with strong resistance to tearing
- **Versatile** – a variety of viscosities satisfy individual demands of clinicians

• 4 - 50ml Cartridge Pack, Instructions, 16 Mixing Tips*

• 32 - 50ml Cartridge Bulk Pack, Instructions, 132 Mixing Tips*

• Putty Catalyst & Base 1 - 230ml Jar each, Instructions, Accessories

Correct Plus Unflavored Impression Cartridges

DESCRIPTION	4 CARTRIDGES	32 CARTRIDGES	PUTTY
Thick'n Thin™ Light Body	REF/PART NO. Q34GA	REF/PART NO. Q34GC	-
Light Body	REF/PART NO. Q34BA	REF/PART NO. Q34BC	-
Universal™ Body	REF/PART NO. Q34DA	REF/PART NO. Q34DC	-
Auto-Mix Putty	REF/PART NO. Q34HA	REF/PART NO. Q34HC	-
Tray Material	REF/PART NO. Q34JA	REF/PART NO. Q34JC	-
Putty - Base (230ml) & Catalyst (230ml)	-	-	REF Q34H

Correct Plus Flavored Impression Cartridges

DESCRIPTION	4 CARTRIDGES	32 CARTRIDGES	PUTTY
Berry Tray Material	REF/PART NO. Q34JAA	REF/PART NO. Q34JAC	-
Berry Auto-Mix Putty	REF/PART NO. Q34HAA	REF/PART NO. Q34HAC	-
Berry Putty - Base (230ml) & Catalyst (230ml)	-	-	REF Q34HCA

IMPRESSION



THE DENTAL
ADVISOR
++++½

Dental PRODUCT
SHOPPER
BEST PRODUCT
2013

Correct Plus™ and Correct Quick™ Bite Registration Material

Correct Plus Bite Superfast and *Correct Quick* Bite are designed to consistently provide highly accurate and detailed bite registrations. *Correct Plus* Bite Superfast and *Correct Quick* Bite have a mousse-like viscosity with a high end hardness/rigidity that facilitates trimming for the perfect bite registration every time. *Correct Plus* Bite Superfast offers a 15-second working time and 1:15-minute total set time, while *Correct Quick* Bite offers a 30-second working time with 1:30-minute total set time.

- 4 - 50ml Cartridge Pack, Instructions, 16 Mixing Tips*
- 32 - 50ml Cartridge Bulk Pack, Instructions, 132 Mixing Tips*

DESCRIPTION	4 CARTRIDGES	32 CARTRIDGES
Correct Plus Bite Superfast Material	REF Q34IA	REF Q34IC
Correct Quick Bite Material	REF Q31E	REF Q31EB

Correct VPS™ Putty


Correct VPS putty is an easy-to-mix, non-oily, heavy viscosity dental putty with a high-end hardness to resist distortion and provide maximum hydraulic drive of the light body materials.

DESCRIPTION	REF/PART NO.
Putty - Base & Catalyst (270ml each)	REF Q03









IMPRESSION


Alginate Material Comparison Chart

Alginate	Description	Indication	Product	Mixing Time	Work Time	Total Set Time	Color
IMAGE	Alginate impression a smooth, glass-like surface that's clear to decipher in the vibrant green color providing smooth detail. The alginate is easy to mix by hand.	Q & FA	Image Dust-Free Regular Set Alginate	0:30	2:20	3:30	
			Dust-Free Fast Set Alginate	0:30	1:45	2:20	

VPS Impression Material Comparison Chart

Viscosity	Description	Indication	Product	Work Time	Intraoral Time	Total Set Time	Color
Thick'N Thin™ Light Body	Thixotropic light-viscosity material that initially stacks on the tooth when syringed and then flows into the sulcus to capture critical subgingival detail.	W	Correct Plus™ Material	1:00	3:30	4:30	
			Correct Plus™ Fast Set Material	0:45	1:45	2:30	
Light Body	A traditional light-viscosity wash material characterized by adequate flowability.	W	Correct Plus™ Material	1:00	3:30	4:30	
Universal™ Body	All-purpose impression material. Use it as a tray material and/or light body. A substitute for alginate for night guards, bleach trays, or any other impression procedure.	W & C	Correct Plus™ Material	1:00	3:30	4:30	
			Correct Plus™ Fast Set Material	0:45	1:45	2:30	
Tray Material	A heavy-viscosity material with a mousse-like consistency. Ideal for Three-in-One impression trays.	Q, FA, TIO	Correct Plus™ Material (Unflavored/Berry Flavored)	1:00	3:30	4:30	
Auto-Mix Putty	A heavy viscosity material in a convenient auto-mix cartridge.	Q, FA, TIO	Correct Plus® Material (Unflavored/Berry Flavored)	1:00	3:30	4:30	
			Correct Plus® Fast Set Material	0:45	1:45	2:30	
Putty	Heavy-viscosity putty	Q & FA	Correct Plus® Material (Unflavored/Berry Flavored)	1:00	3:30	4:30	
			Correct VPS® Material	1:00	3:30	4:30	

Bite Registration Material Comparison Chart

Viscosity	Description	Indication	Product	Work Time	Intraoral Time	Total Set Time	Color
Bite Material/ Bite Superfast	A mousse-like viscosity and high-end hardness. When set, its rigidity facilitates trimming.	Bite Registration	Correct Plus™ Material	0:15	1:00	1:15	
			Correct Quick™ Material	0:30	1:00	1:30	

Indication Abbreviations: W=Wash, C=Custom Tray, Q=Quadrant Tray, FA=Full Arch Tray, TIO=Three-in-One Impression Tray

IMPRESSION

VanR Reversible Hydrocolloid Impression Material



Features & Benefits:

- **Hydrophilic** – works efficiently in the presence of moisture with strong resistance to tearing
- **Versatile** – a variety of viscosities satisfy individual demands of clinicians



Heavy Bodied - Green

The standard by which all impressions are judged and the most popular type. Tempers for 7 minutes at 110°F/43°C and 3 minutes for partial denture impressions. Cool 5 minutes in the mouth, with room temperature water (72°F/23°C). Storage after liquefaction: 5 days.

DESCRIPTION	REF/PART NO.
ECONOPAK (144 Poly Tubes)	11014
12 Boxes of One Dozen Poly Tubes	11011
1 Dozen Tubes (60g ea.)	11012



Qwik Tray Brown

Qwik is a 3 minute gel (40% faster gel). It has easy-to-read colors and a patient pleasing cocoa aroma. Qwik has a tempering time of 3 minutes at 110°F/43°C. Cool 3 minutes in the mouth, with room temperature water (72°F/23°C). Storage after liquefaction: 3 days.

DESCRIPTION	REF/PART NO.
ECONOPAK (144 Poly Tubes)	11403
12 Boxes of One Dozen Poly Tubes	11404
1 Dozen Tubes (60g ea.)	11401



Slate -Pacific Blend

DESCRIPTION	REF/PART NO.
ECONOPAK (144 Poly Tubes)	11302
12 Boxes of One Dozen Poly Tubes	11301
1 Dozen Tubes (60g ea.)	11300



Agarloid - Blue (Fast Tempering)

Tempers for 4 minutes at 110°F/43°C. Cool 5 minutes in mouth, with room temperature water (72°F/23°C). Storage after liquefaction: 5 days.

DESCRIPTION	REF/PART NO.
ECONOPAK (144 Poly Tubes)	11502
12 Boxes of One Dozen Poly Tubes	11503
1 Dozen Tubes (60g ea.)	11501

IMPRESSION



Acculoid - Lavender (Extra Strength)

A heavy bodied formula, that tempers for 5 minutes at 110°F/43°C. Cool 7 minutes in the mouth, with room temperature water (72°F/23°C).

DESCRIPTION	REF/PART NO.
ECONOPAK (144 Poly Tubes)	11111
12 Boxes of One Dozen Poly Tubes	11112
1 Dozen Tubes (60g ea.)	11110



Regular - Red

DESCRIPTION	REF/PART NO.
ECONOPAK (144 Poly Tubes)	11004
12 Boxes of One Dozen Poly Tubes	11002
1 Dozen Tubes (60g ea.)	11001



Backloading Rev. Hydrocolloid Syringe Material

(Box of 1 dozen 25g poly tubes)

Syringe formula Rev. Hydrocolloid packaged in a small poly tube. The nozzle backloads into the 5/16" and the 3/8" Backloading Syringes. Storage after liquefaction: 5 days.

DESCRIPTION	REF/PART NO.
Pink, light viscosity (for displacement use)	11030
White, ultra light (extra working time for multiple preparations)	11031
Blue, heavy viscosity (for placement use)	11034



Qwik - Periwinkle Blue

DESCRIPTION	SHADE	REF/PART NO.
Box of 12 Tubes (25g ea.)	Qwik Wash - Periwinkle Blue	11412
Box of 48 Cartridges (1.8cc ea.)	Syringe Material - Periwinkle Blue	11422



Wash Rev. Hydrocolloid

(Box of 12 tubes, 25g ea.)

One step technique for making Rev. Hydrocolloid impressions; eliminates the need for a syringe for full-coverage crown impressions. Maximum tempering time is 7 minutes. Storage after liquefaction: 5 days.

DESCRIPTION	SHADE	REF/PART NO.
Wash Rev. Hydrocolloid poly tubes		11070
Qwik Wash	Periwinkle Blue	11412

IMPRESSION



CartriLoids, Rev. Hydrocolloid Syringe Material

(1.8cc anesthetic-type glass vials)

Pre-filled syringe material for up to four preparations. To be used with a standard non-aspirating, anesthetic-type syringe. Storage after liquefaction: 5 days.

DESCRIPTION	SHADE	REF/PART NO.
Box of 48	White, Ultra Light (For Displacement Use)	11043
Box of 288	White, Ultra Light (For Displacement Use)	11061
Box of 48	Pink, Light Viscosity (For Displacement Use)	11040
Box of 288	Pink, Light Viscosity (For Displacement Use)	11060
Box of 48	Blue, Heavy Viscosity (For Placement Use)	11038
Box of 288	Blue, Heavy Viscosity (For Placement Use)	11062
Box of 48 (3 day storage)	Periwinkle Blue, fast-gelling (For Use with Qwik Tray Material)	11422



CartriLoids, SuperCartriLoid (Rev. Hydrocolloid)

(Box of 24 large 3.6cc vials)

Pre-filled syringe vials for up to eight preparations. Requires a SuperCartriLoid Syringe. Utilizes the same needle as regular CartriLoid syringes. Storage after liquefaction: 5 days.

DESCRIPTION	SHADE	REF/PART NO.
CartriLoids, SuperCartriLoid	Pink, Light Viscosity (For Displacement Use)	11065
CartriLoids, SuperCartriLoid	Blue, Heavy Viscosity (For Placement Use)	11068

IMPRESSION



Syringe Sticks

Rev. Hydrocolloid Syringe Material

Prefabricated Syringe Material which is placed in either the 5/16" or 3/8" diameter backloading Rev. Hydrocolloid syringes. Storage after liquefaction: 5 days

Syringe Sticks		Rev. Hydrocolloid Syringe Material
DESCRIPTION	SHADE	REF/PART NO.
5/16" Diameter Hydrocolloid Sticks	Pink, Light Viscosity (For Displacement Use)	11016
3/8" Diameter Hydrocolloid Sticks	Pink, Light Viscosity (For Displacement Use)	11018
3/8" Diameter Hydrocolloid Sticks	Blue, Heavy Viscosity (For Placement Use)	11098

CartriLoid/SuperCartriLoid Basket		Holds 8 CartriLoids or 5 SuperCartriLoids
DESCRIPTION		REF/PART NO.
Metal CartriLoid Basket		11058

CartriLoid Needles, Blunt-end, 18 Gauge Disposable		Box of 100 needles for CartriLoids and SuperCartriLoids
DESCRIPTION	SHADE	REF/PART NO.
18 gauge x 3/8" (100 per package)	White	11049

CartriLoid and SuperCartriLoid Syringes		(breech loading, non-aspirating)
DESCRIPTION		REF/PART NO.
Metal CartriLoid Syringe		11053
Metal SuperCartriLoid Syringe		11067

Syringe, Backloading Rev. Hydrocolloid		18 gauge backloading
DESCRIPTION		REF/PART NO.
Accommodates 5/16" Stick Syringe or Wash Rev. Hydrocolloid		12018
Accommodates 3/8" Stick Syringe or Wash Rev. Hydrocolloid		12041

IMPRESSION



Extend-A-Pour

An impression storage solution designed for use with Rev. Hydrocolloid and alginate. Enables the office to delay impression pouring for up to 4 weeks after the impression is immersed.

DESCRIPTION	REF
12oz Bottle	13005



B-Tron Hydroprocessor, 110 Volt

Modular design with three interchangeable plug-in units which contain all of the electronic components, microcomputers and heating elements for easy, convenient servicing.

Separate display panels indicate the temperature of each bath, while function lights indicate Liquefying and Timing modes. There are also special displays to warn of power interruptions and when water levels become too low. The *B-Tron* is factory programmed with preset temperatures and is ready to use.

Dimensions (B-Tron 301):

- **Width** – 16 1/4"
- **Depth** – 10" (plus, allow 1 1/2" for power supply cord)
- **Height** – 10 1/2" (lids closed), 16" (lids open)
- **Weight** – 15 lb
- **Electrical** – 120v 50-60 Hz, 10.4 amps

Also available in foreign operating voltages

B-Tron 301, 110v (60 Hz)

DESCRIPTION	REF/PART NO.
B-Tron 301, 110v (60Hz) (3 bath unit: Boil/Auxiliary Store, Store and Temper)	15510

B-Tron Parts

DESCRIPTION	REF/PART NO.
B-Tron Plug-in Electronic Module (110v)	15550
B-Tron/Accutron Cosmetic Refurbish Kit (includes 3 tanks, 3 tank lids, 3 baskets and 3 control panel labels)	15590

IMPRESSION



Tacky Stops

Disposable, self-adhesive, occlusal stops (1/16" x 3/8" x 1/16") which adhere to the occlusal surface of trays. Tacky Stops prevent over-seating of tray.

DESCRIPTION	REF/PART NO.
Box of 300	12051



Post Damming Plastic

Used to contain Reversible Hydrocolloid or any impression material in a maxillary tray and to increase compression for maximum displacement during seating of the tray. Tray must be clean and dry with no trace of glove powder for best adherence to the tray.

DESCRIPTION	REF/PART NO.
3/4 lb box (56 pre-cut strips, 3" x 1/2" x 1/8")	13061



Periphery Plastic

Designed to augment tray height. Also used along the periphery edge of trays for patient comfort or for a peripheral roll. Tray must be clean and dry with no trace of glove powder for best adherence to the tray.

DESCRIPTION	REF/PART NO.
3/4 lb box (49 pre-cut strips, 5" x 3/4" x 1/8")	13060



PrepWet

A lubricant which is sprayed over the entire arch after the retraction cord is removed. PrepWet allows the Reversible Hydrocolloid syringe material to flow freely throughout the sulcus, and easily remove the impression tray. PrepWet reduces tearing when sprayed on mouth, prior to seating the alginate impression.

DESCRIPTION	REF/PART NO.
PrepWet Kit (1 pkg. of 3, 4 oz sprayer bottles w/100 Dispenser Straws)	13097
Dispenser Straw (pkg. of 100)	13096

CEMENTS

Mojo™ Veneer Cement

Mojo Veneer Cement is a light cure, esthetic cementation system that is designed for use with porcelain, ceramic, and composite veneer restorations. The simple, highly versatile selection of shades allows this material to be used with a wide range of veneer cases, while offering the least amount of shift available.* The corresponding water-soluble Try-In gels, included in the comprehensive kit, allow for a perfect match every time.

Features and Benefits

- **Easy cleanup** – simply light cure for 1–2 seconds and remove excess cement
- **Undetectable shade shift** – the restoration stays the same shade as on the first day you cured the cement
- **Versatile shade options** – allows you to warm or brighten the restoration easily for a natural look; shades can be mixed for superlative shade customization
- **Simple storage** – no refrigeration needed and easily stored at room temperature



REF N06K **Mojo Veneer Cement Kit**

- 4 - Syringes: Cement (2.3gm) Clear, Dark, Light, *Movie Star White*
- 4 - Syringes: Try-In (1.6gm) Clear, Dark, Light, *Movie Star White*
- 1 - 4ml Bottle *Bond-1* Primer/Adhesive
- 1 - 1ml Syringe 37% Phosphoric Acid Etching Gel
- 1 - 3ml Bottle Silane
- 40 - Elephant Tips
- Instructions

Mojo Veneer Cement Refills

Includes 1 - 2.3gm Syringe, 5 - Elephant Tips, Instructions

DESCRIPTION	REF/PART NO.	DESCRIPTION	REF/PART NO.
<i>Mojo Light</i>	N06A	<i>Mojo Clear</i>	N06C
<i>Mojo Dark</i>	N06B	<i>Mojo Movie Star White</i>	N06D

Mojo Try-In Gel Refills

Includes 1 - 1.6gm Syringe, 5 - Elephant Tips, Instructions

DESCRIPTION	REF/PART NO.	DESCRIPTION	REF/PART NO.
<i>Mojo Light</i>	N06AA	<i>Mojo Clear</i>	N06AC
<i>Mojo Dark</i>	N06AB	<i>Mojo Movie Star White</i>	N06AD

Mojo Accessories

DESCRIPTION	REF/PART NO.
100 - Elephant Tips	N06E

*Compared to leading materials.

CEMENTS



Breeze™ Self-Adhesive Resin Cement

Breeze Self-Adhesive Resin Cement is specifically designed to make the cementation of crowns, bridges, inlays, onlays and posts faster and easier by eliminating individual etching, priming, bonding, and mixing steps. Formulated with advanced resin technology, *Breeze* Self-Adhesive Resin Cement provides one of the strongest retention available in a self-adhesive cement with the quick and easy to use auto-mix delivery system. Upon application, *Breeze* Cement quickly goes to work to condition dentin, enamel, and the restoration all in one quick and simple step. By eliminating the etching, priming and bonding steps, both postoperative sensitivity and procedure time are drastically reduced improving peace of mind and freeing up precious chair time.

Features & Benefits:

- **Exceptional bond strength** - adheres to dentin and enamel without a bonding agent
- **Convenient auto-mix delivery system** - single-step accurate placement with minimal waste
- **Versatile shade selection** - choose from A2, Translucent, and Opaque White to satisfy all cosmetic needs
- **Dual Cure** - light cure or self cure for maximum versatility

REF N97 *Breeze* Self-Adhesive Resin Cement Intro Kit

3 - 4ml Syringes A2	40 - Auto-mix Tips
1 - 4ml Syringe Translucent	10 - Intra-oral/Auto-mix Tips
1 - 4ml Syringe Opaque White	• Instructions
1 - 3ml Bottle Silane	

Breeze Cement Refills

Includes 1 - 4ml Syringe, 12 Mixing Tips, Instructions

DESCRIPTION	REF/PART NO.
A2	N97B
Translucent	N97A
Opaque White	N97E

Accessories

DESCRIPTION	REF/PART NO.	DESCRIPTION	REF/PART NO.
1 - 3ml Bottle Silane	N04EA	50 - Intra-oral/Auto-mix Tips	N72AA
50 - Auto-mix Tips	N33CA	100 - Intra-oral/Auto-mix Tips	N72AB
100 - Auto-mix Tips	N33C	1 - 1:1 Auto-Mix Syringe Dispenser Gun	N75S

CEMENTS

Cement-It™ Universal Crown & Bridge Resin Cement

Cement-It Universal Crown & Bridge Material is a dual cure resin cement for all of your indirect bonding needs. The system delivers high bond strengths, reduced sensitivity, and versatile bonding for any restoration. *Cement-It* is available in a quick, easy-to-use auto-mix delivery system and is available in three shades to satisfy virtually all cosmetic needs: A2, Opaceous White, and Translucent.

Features & Benefits:

- **Convenient** – auto-mix syringe eliminates mixing guesswork and delivers consistent results
- **Versatile** – compatible with all crowns, bridges, inlays, onlays, and posts
- **Durable** – creates long lasting restorations with bond strengths up to 25 MPa



Cement-It Universal Resin Cement C&B Refills

Includes 1 - 4ml Syringe, 12 Mixing Tips, Instructions

DESCRIPTION	REF/PART NO.
1 - 4ml Syringe A2 Catalyst & Base	N33A
1 - 4ml Syringe Op. White Catalyst & Base	N33AB
1 - 4ml Syringe Translucent Catalyst & Base	N33AC

Accessories

DESCRIPTION	REF/PART NO.	DESCRIPTION	REF/PART NO.
1 - 3ml Bottle Silane	N04EA	50 - Intra-oral/Auto-mix Tips	N72AA
50 - Auto-mix Tips	N33CA	100 - Intra-oral/Auto-mix Tips	N72AB
100 - Auto-mix Tips	N33C	1 - 1:1 Auto-Mix Syringe Dispenser Gun	N75S

BONDING AGENTS



Bond-1™ Primer/Adhesive

The *Bond-1* Primer/Adhesive System is a 5th generation total-etch bonding agent that provides high bond strength in a convenient single-bottle application. The versatile system bonds both direct and indirect restorations. *Bond-1* creates a sealed surface with an intricate network of polymerized primers and unfilled resin. This results in a high bond strength to dentin of up to 31.0 MPa and a film thickness of only 8µm. Its ease of application—combined with its low film thickness—provides the ideal surface to bond to all types of restorative materials.

Features & Benefits:

- **Flexible** – bonds to all direct and indirect restorations
- **Easy-to-use** – total-etch bonding technique that's familiar and proven
- **Resilient** – bond strengths up to 31 MPa to dentin / 27 MPa to enamel

REF N011 *Bond-1* Primer/Adhesive Kit

- 1 - 4ml Bottle *Bond-1* Primer/Adhesive
- 1 - 3ml Bottle *Bond-1* Dual Cure Activator
- 2 - 5ml Syringes 37% Phosphoric Acid Etching Gel
- Instructions, Accessories

Bond-1 Refills

Includes Instructions, Accessories

DESCRIPTION	REF/PART NO.
1 - 4ml Bottle <i>Bond-1</i> Primer/Adhesive	N011AA
1 - 6ml Bottle <i>Bond-1</i> Primer/Adhesive	N011AB
1 - 3ml Bottle <i>Bond-1</i> Dual Cure Activator	N011AC
2 - 5ml Syringes 37% Phosphoric Acid Etching Gel	N011B

BONDING AGENTS



Bond-1™ SF Solvent-Free SE Adhesive

Bond 1-SF is a self-etch, one-coat bonding agent that is solvent-free — eliminating the need to air dry. This feature protects against sensitivity and prevents evaporation for stability and consistency. *Bond-1 SF* is designed for direct composite restorations and dual cure materials.* *Bond-1 SF* provides optimal bond strengths up to 30.4 MPa in just three basic steps: apply evenly, rub for 20 seconds, and light cure.

Features & Benefits:

- **Solvent-free** – protects against sensitivity, eliminating common technique-sensitive issues such as under- and over-drying; also prevents evaporation
- **No air drying** – saves you precious chair time; doesn't leave you guessing "how dry is too dry"
- **Self-etch** – no need to acid etch; saves time and protects against sensitivity
- **Versatile** – can be used in conjunction with light cure and dual cure materials*
- **Unique handling properties** – enables an even spread of the material, ensuring complete coverage and optimal results
- **One coat** – easy to use; reduces procedure time

Bond-1 SF Solvent-Free SE Adhesive Syringe Refills

DESCRIPTION	REF/PART NO.
2 - 1ml Syringes, 20 - 25 Gauge Flocked Needle Tips, Instructions	N03N

Bond-1 SF Accessories

DESCRIPTION	REF/PART NO.
100 - 25 Gauge Flocked Needle Tips	N03T
100 - Disposable Applicators	N10B

*Compatible with dual cure materials only when the initial layer of dual cure material is light cured.

BONDING AGENTS



Bond-It™ Enamel / Dentin Bonding System

Bond-It is a 4th generation total-etch bonding system that offers a reliable technique proven over time. It is a perfect choice for bonding dual and self cure composites as well as amalgam. Known as the “golden standard” of bonding, this familiar technique assures the ideal surface to bond any restorative material, including direct and indirect composite, porcelain, metal, and amalgam.

Features & Benefits:

- **Flexible** – bonds to all surfaces
- **Familiar** – easy to use, reliable system
- **Ideal** – perfect for bonding under amalgam

REF N01 *Bond-It Bonding System Kit*

- 1 - 6ml Bottle Primer A
- 1 - 6ml Bottle Primer B
- 2 - 3ml Syringe 37% Phosphoric Acid Etching Gel
- 1 - 6ml Bottle Resin Light Cure
- 1 - 3ml Bottle Resin Dual Cure Activator
- 20 - 22 Gauge Tips
- Instructions, Accessories




Bond-It Enamel / Dentin Bonding Refills

Includes Instructions, Accessories

DESCRIPTION	REF/PART NO.
1 - 6ml Bottle Primer A	N01A
1 - 6ml Bottle Primer B	N01B
1 - 6ml Bottle each Primer A and B	N01C
1 - 6ml Bottle Resin Light Cure	N01D
1 - 3ml Bottle Resin Dual Cure Activator	N01E
1 - 6ml Bottle Resin Light Cure and 1 - 3ml Bottle Dual Cure Activator	N01F
2 - 3ml Syringes 37% Phosphoric Acid Etching Gel, 10 - 22 Gauge Tips	N01H

ACCESSORIES

Etching Gel




PRODUCT PHOTO	DESCRIPTION	QUANTITY	REF/PART NO.
	37% Phosphoric Acid Etching Gel – 1ml/1.3gm Syringe 20 - 22 Gauge Tips	2	T06
	37% Phosphoric Acid Etching Gel – 3ml/4gm Syringe 20 - 22 Gauge Tips	2	N01H
	37% Phosphoric Acid Etching Gel – 5ml/6.9gm Syringe 20 - 22 Gauge Tips	2	N01IB

Cement/Adhesive Accessories

Adhesives Accessories	
DESCRIPTION	REF/PART NO.
1 - 3ml Bottle Silane	N04EA
100 - Disposable Applicators	N10B
1 - 1:1 Auto-Mix Syringe Dispenser Gun	N75S
100 - 25 Gauge Flocked Needle Tips	N03T

Impression Accessories

Impression Accessories			
DESCRIPTION	REF/PART NO.	DESCRIPTION	REF/PART NO.
1 - 1oz. Bottle Impression Tray Adhesive	Q12	1 - Impression Dispenser Gun	Q19
2 - 1oz. Bottles Impression Tray Adhesive	Q12A		

PRODUCT PHOTO	DESCRIPTION	QUANTITY	REF/PART NO.
	37% Phosphoric Acid Etching Gel – 1ml/1.3gm Syringe 20 - 22 Gauge Tips	2	T06
	37% Phosphoric Acid Etching Gel – 3ml/4gm Syringe 20 - 22 Gauge Tips	2	N01H
	37% Phosphoric Acid Etching Gel – 5ml/6.9gm Syringe 20 - 22 Gauge Tips	2	N01IB


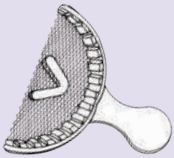
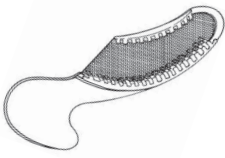

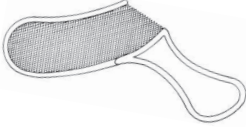
ACCESSORIES

Impression Trays

Clear Impression Trays (Excellent to use with TempSpan™ Clear Matrix)

	DESCRIPTION	QUANTITY	REF/PART NO.
	Upper Right / Lower Left Clear Trays	12	Q17F
	Upper Left / Lower Right Clear Trays	12	Q17G
	Anterior Clear Trays	12	Q17H

Three-In-One Impression Material Trays

	Full Arch Tray	30	Q17A
	Anterior Tray	35	Q17B
	Quadrant Tray	35	Q17C
	Posterior Tray	50	Q17D
	Sideless Posterior Tray	50	Q17E


ACCESSORIES

Tip Refills

		DESCRIPTION	QUANTITY	REF/PART NO.
Bond-1™ SF Solvent-Free SE Adhesive		25 Gauge Flocked Needle Tip	100	N03T
Fusio™ Self-Adhesive Flowable Composite		19 Gauge Bent Needle Tip	50	N28NT
Flow-It™ ALC Flowable Composite		20 Gauge Bent Needle Tip	50	T26A
			100	T26C
Phosphoric Acid Etching Gel		22 Gauge Bent Needle Tip	100	T26B
Build-It™ FR Core Material		Wastesaver Tips (Blue)	48	N32E
Build-It™ FR Core Material		Intra-oral Tips (Yellow)	48	N32D
Build-It™ FR Mini-Mix Core Material		Mixing Tips & Intra-oral Tips	48	N32GA
			100	N32GB
TempSpan™ Temporary C&B Material		Mixing Tips	30	N69AH
			100	N69AI
ZONE, ZONEfree, Embonte Temporary Cements		Automix Syringe Tips	8	27048

ACCESSORIES

Tip Refills


		DESCRIPTION	QUANTITY	REF/PART NO.
TempSpan™ Transparent Temporary Cement Breeze™ Self-Adhesive Resin Cement Cement-It™ Universal C&B Resin Cement		Mixing Tips	50	N33CA
			100	N33C
Breeze™ Self-Adhesive Resin Cement Cement-It™ Universal C&B Resin Cement		Intra-oral & Auto-mix	50	N72AA
			100	N72AB
Mojo™ Veneer Cement		Elephant Tips	100	N06E
Impression Materials		Intra-oral Mixing Tips <i>(used with Pink or Green mixing Tips)</i>	100	Q16
Correct Plus™ Impression Materials Light Body, Thick-n-Thin		Pink Mixing Tips	100	Q19AB
TempSpan™ Clear Matrix Impression Materials Bite Registration, Tray Material, Universal Body		Green Mixing Tips	100	Q19BB
Impression Materials Auto-Mix Putty		Purple Mixing Tips	100	Q19CB
Correct Plus™ Mighty Mix™ Vinyl Polysiloxane Impression Material		Dynamic Mixing Tips	50	Q34MM
		Bayonet Rings	2	Q34MR

ACCESSORIES

Dispenser Guns

	DESCRIPTION	REF/PART NO.
	Single-Dose Dispenser Gun <i>Product Lines: Simile, Artiste Nano, Alert, Build-It Light Cure</i>	N75A
	1:1 Auto-Mix Syringe Dispenser Gun <i>Product Lines: Build-It FR Mini-Mix, Breeze, Cement-It, TempSpan Transparent Temporary Cement</i>	N75S
	25ml 1:1 Dispenser Gun <i>Product Lines: Build-It FR</i>	N32B
	50ml 1:1 Dispenser Gun <i>Product Lines: Correct Plus, Correct Quick, TempSpan C&B, TempSpan Clear Matrix Material</i>	Q19

Non-Slip Mixing Pads

	DESCRIPTION	QUANTITY	REF/PART NO.
	MIXING PADS #0-X SMALL	12	31092
	MIXING PADS #1-SMALL	8	31093
	MIXING PADS #2-MEDIUM	4	31094
	MIXING PADS #3-LARGE	2	31095
	MIXING PADS X-LARGE	1	31096

RETURNED MERCHANDISE

Pentron will honor a 30-day money back guarantee on all products. Please contact your distributor for additional return policies and procedures.

PRODUCT PERFORMANCE COMPLAINTS

In the event of a suspected product performance complaint, contact your dealer or Pentron distributor for a return merchandise authorization (RMA). DEALERS NOTE: Pentron will not issue credit or product replacement unless all required complaint information is received.

LIMITED WARRANTY- LIMITATION OF PENTRON'S LIABILITY

Pentron's technical advice, whether verbal or in writing, is designed to assist dentists in using Pentron product. The dentist assumes all risk and liability for damages arising out of improper use of Pentron product. In the event of a defect in material or workmanship, Pentron's liability is limited, at Pentron's option, to replacement of the defective product or part thereof, or reimbursement of the actual cost of the defective product. In order to take advantage of this limited warranty, the defective product must be returned to Pentron. In no event shall Pentron be liable for any indirect, incidental, or consequential damages.



PENTRON

Toll Free: 800.551.0283 | Phone: 714.516.7557 | Fax: 877.677.8844
www.pentron.com | www.pentronstrada.com

Pentron Catalog_Rev02